



PHILIPS

RESPIRONICS


SimplyGo Mini

Portable oxygen
concentrator

Reliability that's
simply portable

PHILIPS
RESPIRONICS

Built for the demands of today's
active patients

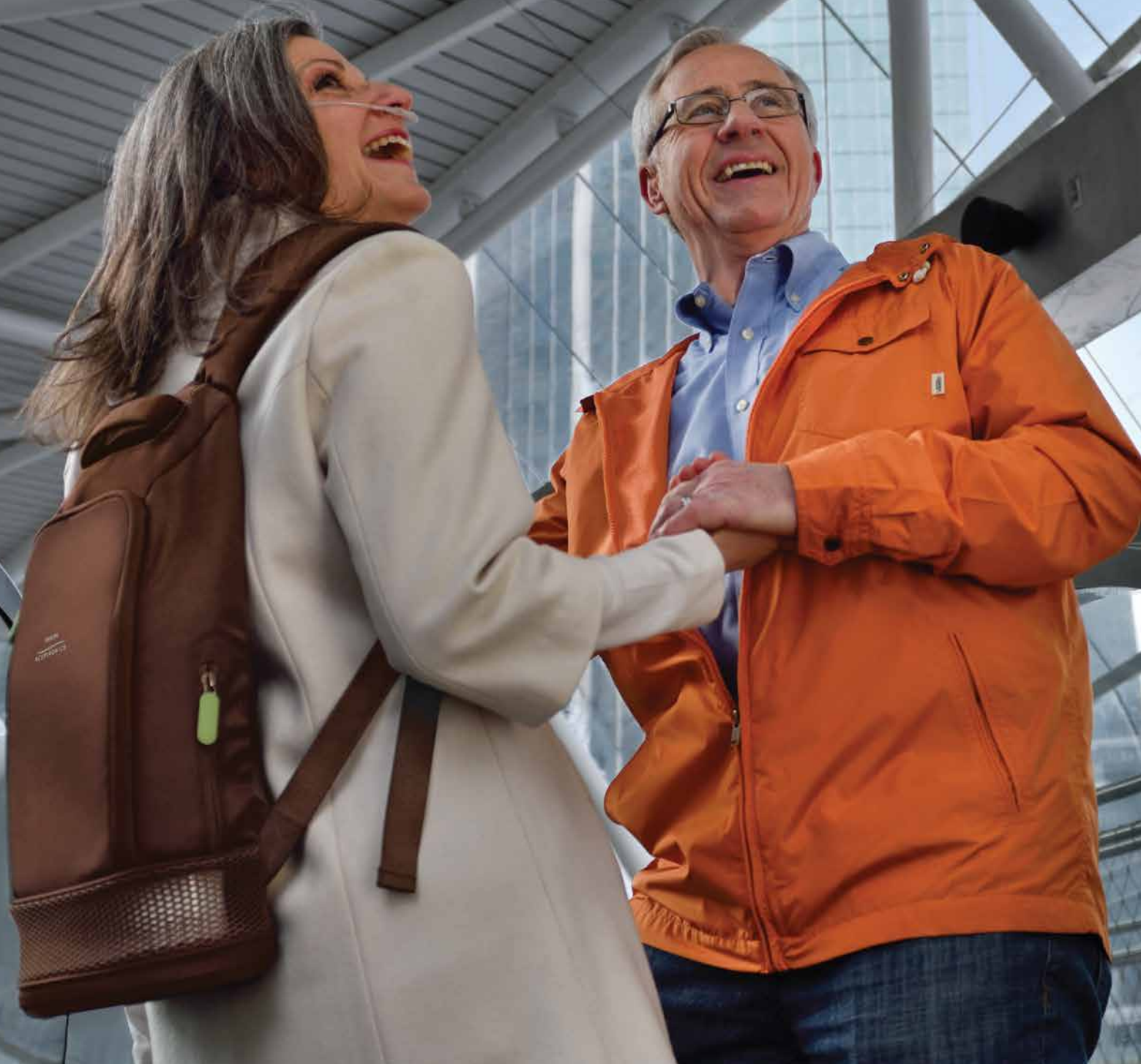

The SimplyGo Mini portable oxygen concentrator (POC) is designed to help you grow your business by giving you the ability to offer light weight, attractive features, and true reliability that today's patients demand in a POC.

At Philips Respironics, we have earned our reputation for designing and building excellence into all of our products. Now, with SimplyGo Mini, you can offer your patients the choice for excellence in lightweight portable oxygen concentrators.

Reliability matters to your business

When you place a POC with a patient, you need to know that it's going to have the reliability and ruggedness to help your business:

- Eliminate recurring oxygen delivery costs and cylinder hassles
- Reduce maintenance costs
- Ease inventory management
- Minimize patient calls

SimplyGo Mini has undergone the same rigorous testing regimen as that of our SimplyGo POC. SimplyGo Mini was subjected to real-world conditions including impacts, vibrations, and temperatures to deliver reliable performance day in and day out.

Simply smart design

As the smallest and lightest POC ever developed by Philips Respironics, SimplyGo Mini is built for today's patients with features like:

- Attractive and sleek design
- Easy-to-remove external battery
- Light weight (5.0 lbs.)
- Intuitive, easy-to-read screen
- Strong, comfortable carrying case



The manufacturer of this POC has determined this device conforms to all applicable FAA requirements for POC carriage and use on board aircraft.



Complies with
RTCA/DO-160G
section 21, category M

Why purchasing a quality product matters



Potential increase in patient retention



Potential increase in clinician confidence and referrals



Increase in employee value-added activities because of fewer failure calls



Decrease in inventory because of fewer swap-outs



Decrease in resources and overtime payments to deal with failed equipment calls and pick-ups



Decrease in day-to-day business disruptions because of product failures

One COPD patient care solutions

An integrated care model that puts the patient first
One COPD takes a patient-centric approach—based on clinical research and backed by Philips devices, support, and expertise.

Technology, education, and connectivity all work together to provide a seamless care experience geared towards patient satisfaction, better patient care, and reduced cost of care.





Product specifications

Oxygen concentrations ¹	At least 87% at all settings (maximum of 96%) ¹ Based on atmospheric pressure of 14.7 psia (101 kPa) at 70° F (21° C)		
Outlet pressure	20 psig		
Flow settings and pulse volumes	Setting	Nom. pulse vol. at 20 BPM	Max. minute vol. output
	1	11 ml	220 ml/min
	2	22 ml	440 ml/min
	3	33 ml	660 ml/min
	4	44 ml	880 ml/min
	5	50 ml	1,000 ml/min
	+/- 25% or 6 ml, whichever is greater (average of 20 consecutive pulses) over the rated environmental range		
Sound level	43 dBA typical at settings 2 and 20 BPM; 46 dBA typical at settings 5 and 20 BPM		

Physical characteristics

Size	With standard battery: 9.4 in. x 8.3 in. x 3.6 in. (23.9 cm x 21.1 cm x 9.1 cm) With extended battery: 10.2 in x 8.3 in x 2.6 in (25.9 cm x 21.1 cm x 9.1 cm)		
Weight	5 lbs. (2.3 kg) with standard battery installed 6 lbs. (2.7 kg) with extended battery installed		
User interface	Touch screen LCD		

Operating conditions

Operating temperature	41° F to 95° F (5° C to 35° C)		
Relative humidity	15% to 93%		
Altitude	Up to 10,000 ft. (3,048 m)		

AC power

Input	100 to 240 VAC, 50/60 Hz		
Output	19 VDC, 6.3 A max		
Power consumption	120 W while charging		

DC power

Input	12 – 16 VDC		
Output	19 +/- 5% VDC, 6.3 A		

Audible alarm

Sound level	56 dBA or greater when measured at 1 m		
Alarm type	No breath, low oxygen concentration, high breath rate, technical fault, low battery, no flow, depleted battery, external power failure		

Battery

Chemistry	Lithium ion 14.4 VDC (nominal)		
Battery options	Standard battery, Extended battery (equals two standard batteries)		
Dimensions	Standard battery: 7.9 in. x 3.6 in. x 1.0 in. (20.1 cm x 9.1 cm x 2.5 cm) Extended battery: 7.9 in. x 3.6 in. x 1.9 in. (20.1 cm x 9.1 cm x 4.8 cm)		
Weight	Standard battery: 1.1 lbs. (0.50 kg) Extended battery 2.1 lbs. (0.95 kg)		
Watt-hour rating per standard battery	3.4 Ah x 3.6V x 8 cells = 97.9 Wh per standard battery		
Charge time	Standard battery: 4 hours maximum Extended battery: 8 hours maximum		

Typical battery duration on pulse setting 2, 20 BPM²

Standard battery	Up to 4.5 hours		
Extended battery	Up to 9 hours		

²Battery duration times are based on new, fully charged batteries used with a new SimplyGo Mini System. Battery duration times will degrade with battery age, environmental operating conditions, use over time, and operational condition of the concentrator.

Ordering information

Device	Part number
SimplyGo Mini system with standard battery includes: SimplyGo Mini portable oxygen concentrator, carrying case, one standard rechargeable lithium ion battery, DC and AC power cords, accessory bag, and user manual	US: 1113601 International: 1113604
SimplyGo Mini system with extended battery includes: SimplyGo Mini portable oxygen concentrator, carrying case, one extended rechargeable lithium ion battery, DC and AC power cords, accessory bag, and user manual	US: 1113602 International: 1113605



Replacement accessories	Part number
Standard lithium ion battery	1116816
Extended lithium ion battery	1116817
SimplyGo Mini AC power supply	1116818
SimplyGo Mini AC power cord	1082664
SimplyGo Mini DC power supply	1116819
SimplyGo Mini external battery charger	1116830
User Manual	1118724

Accessories - Bags - Brown	Part number
SimplyGo Mini carry bag and strap, brown	1119897
SimplyGo Mini accessory bag, brown	1116824
SimplyGo Mini backpack, brown	1116836

Accessories - Bags - Black	Part number
SimplyGo Mini carry bag and strap, black	1119928
SimplyGo Mini accessory bag, black	1119893
SimplyGo Mini backpack, black	1119894

Warranties	Part number
SimplyGo Mini 3-year system, 1-year sieve warranty	1120121
SimplyGo Mini 3-year system, 3-year sieve warranty	1122070
SimplyGo Mini 2-year sieve warranty extension ³	1120123

³2-year sieve warranty extension is a post-sale warranty purchase option for those units purchased initially with the 3-year system, 1-year sieve warranty.

Learn more

For more information about SimplyGo Mini, call customer service at 1-800-345-6443.

Visit www.philips.us/simplygomini

The Right Fit

SimplyGo Mini is part of our Right Fit, a complete portfolio of oxygen products and programs inspired by patients and built for business.



© 2016 Koninklijke Philips N.V. All rights reserved.
Specifications are subject to change without notice.

www.philips.com/respironics

Caution: US federal law restricts these devices to
sale by or on the order of a physician.

RRDPGH KW 6/7/16 MCI 4107375 PN 1121655
1010 Murry Ridge Lane, Murrysville, PA 15668
800-345-6443 · 724-387-4000